

Production Engineering - Medical Equipment Division

[PEMED HOME](#) 6035 E. 38th Avenue, Denver, Colorado 80207 USA (303) 393-7800 FAX (303) 393-1482

Pharmacy

To get this entire section as a printable PDF document [click here](#)

Automated tablet dispensing and packaging system - \$19,750



List: \$115,800. MedCon: EX. JV350SL automatic tablet distributing and packing system (also called Autopack), manufactured by JVM Co.. 120 V, 60 Hz.

Features:

- dispenses and packs any combination of up to 350 different capsules
- 350 dispensing canisters arranged in 5 modules
- infrared beam tablet detection system
- thermal transfer package printing system
- self contained packing and printing, from one-pill to multi-pill tablet packing
- printer will mark each package up to 30 lines of information, i.e., drug name, bar code, efficacy, patient name, dosage, etc
- Font, text size and number of letters are adjustable.
- does both unit dose (single pill packing) and multi dose (packing of two or more pills) packages
- 10 inch diagonal touch screen monitor
- computerized system with Microsoft windows 2000 Professional Edition
- circuit board access doors
- special tablet system trays, two included

Removable loading hoppers



Form, fill and seal packaging unit, with loading collection hopper and print mechanism (produces paper or foil packages)



Views of the storage modules from front, right and left sides



Interior layout, and power supply and circuit boards



Manufacturer's label



spare parts package: rotary encoder, circuit board EH0007, Innotech Systems Spitfire with AC adapter, and two rolls of barcode tape



Overall dimensions of the unit:

- width: 43 in. (49 with monitor screen mounted)
- depth: 47.5
- height: 76.5

Shipping weight is about 2124 pounds.

Automatic pipetter/diluter - \$875



List: \$3,925. MedCon: VG. Gilford model 6065 automatic pipettor/diluter. 115 V, 60 Hz.

Features:

- pipetter handle with 2.5 inch shaft
- twin syringes
- continuous/auto/single toggle
- cycle light

Overall dimensions of the unit:

- width: 7 in.
- depth: 6.75
- height: 10.25

Shipping weight is about 10 pounds.

Pharmacy dispensing pump with stand - \$2,500



List: \$9,300. MedCon: VG. Baxa catalog no. 095 peristaltic repeater pump. 120 V, 60 Hz.

Features:

- rolling stand with 24 x 24 inch base and 5 casters (two locking)
- LCD window shows volume, preset, filled and rate
- rate adjusts to 99
- 1 in. rollers
- 1 cm tubing channel
- start/stop button
- reverse button
- plastic lid with flow direction arrows
- power on/off pushbutton with green indicator light
- footswitch included

Overall dimensions of the pump:

- width: 7 in.
- depth: 9.5

- height: 10; on stand the height is 42.5

Shipping weight is about 55 pounds.

Pill counter - \$3,475



List: \$8,700. MedCon: G. C. E. King electronic tablet counter, model TB4 . 120 V, 60 Hz.

Features:

- control panel has batch counter and batch totalizer; counts from 1 to 10,000
- adjustable dials for turntable speed and vibrator
- dust and chip removal
- hopper with lid
- two dispensing ports
- six plastic dispensing spouts (two each large, medium and small openings) are included
- 5.5 foot power cord with hospital grade plug



Overall dimensions of the unit:

- width: 35.5 in.
- depth: 27
- height: 26

Shipping weight is about 300 pounds.

Pill counter - \$975



List: \$2,300. MedCon: G. Kirby Lester model KL-X microprocessor controlled tablet counter. 110 V.

A simple, manually operated pill counters. Tablets and pills are dumped in top of unit, and are counted and passed into a container below.

Features:

- top loading chimney
- test and zero functions
- power on/off rocker switch

Overall dimensions of the unit:

- width: 10.75 in.
- depth: 11.5
- height: 19

Shipping weight is about 20 pounds. 2 units in stock.

Pill counter - \$1,375



List: \$6,700. MedCon: G. C. E. King medication dispenser. 115 V, 60 Hz. Needs hopper assembly and one dispensing chute.

Features:

- control panel has batch counter and batch totalizer
- adjustable dials for turntable speed and vibrator
- two dispensing ports

Overall dimensions of the unit:

- width: 35 in.
- depth: 25
- height: 14

Shipping weight is about 200 pounds.

Automated medication dispensing system - \$1,475



List: \$180,000. MedCon: inc. Baxter ATC Profile automated medication dispensing system. 117 V, 60 Hz. Probably needs refurbishing.


We are selling this unit for parts.

Features:

- 240 medication bins
- dispensing and packaging system
- LCD window shows operation status and advanced diagnostics
- two RS-232C ports



control panel showing various operation features

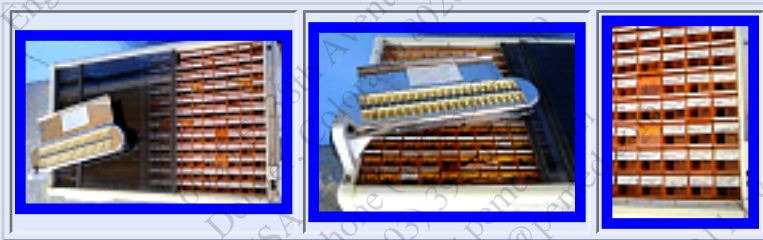

views of the medication bins

Overall dimensions of the unit:

- width: 46 in.
- depth: 29.5
- height: 45

Shipping weight is about 450 pounds.


Unit dose packaging machine sealer - \$1,975



List: \$6,230. MedCon: VG. Euclid Medical Systems The Cadet (X-temp series) strip packaging machine. 115 V. Includes Plexiglas disc for pill feeding and sorting. Does not have the stencil stamping attachment.

Features:

- seals pills between poly-foil and poly-cello
- processes 60 units/min., 1.5 x 2 in.
- 1 cycle and continuous run modes
- power on/off and heater on/off switches with indicator lights
- instruction manual
- supplies include 9 rolls each of poly-foil and poly-cello

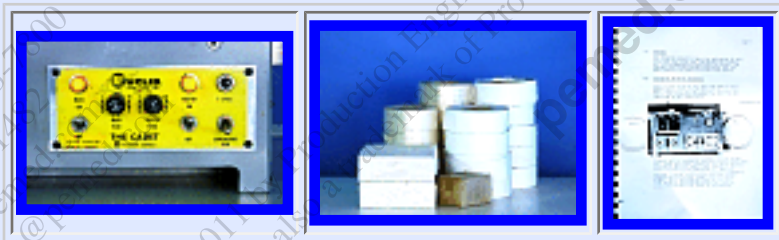

Overall dimensions of the unit:

- width: 22 in.
- depth: 11
- height: 10

Shipping weight is about 70 pounds.

Screw capper and crimper - \$3,750








List: \$16,425. MedCon: EX. Cozzoli Machine Company model SC50T. 115 V, 60 Hz. 1/6 HP GE motor. Provides uniform tightening of screw caps or crimping of aluminum seals. Supports closure sizes of up to 95mm and a production rate of up to 50 per minute by manual placement of containers.

Includes small air/vacuum pump as shown.

Overall dimensions of the unit:

- length: 18 in.
- width: 14
- height: 30

detail of rotary head	
side view	
back view	
Welch/Thomas compressor model 2545B-01	
manufacturer's plate	

Shipping weight is about 130 pounds.

Tube sealer - \$475



List: \$1,925. MedCon: VG. Terumo tube sealer model ACS-152. Seals PVC and other plastic tubing of less than 6mm outside diameter using dielectric heating.

Overall dimensions of unit:

- length: 7 in.
- width: 12
- height: 6

Shipping weight is about 25 pounds.

Packaging/sealing machine - \$975

List: \$2800. MedCon: VG. Doughboy tabletop or countertop package ceiling machine. For ceiling packages and pouches, containing either pills or sterile surgical instruments.




Shipping weight is about 75 lbs.

Tablet packager - \$975



List: \$5,600. MedCon: G. Pharma-Pak tablet packaging machine, model XL. Packages tablets in clear and opaque white tape.

With tray and wheel for small volume packaging. Unit handles a 2" web width. on stand as shown

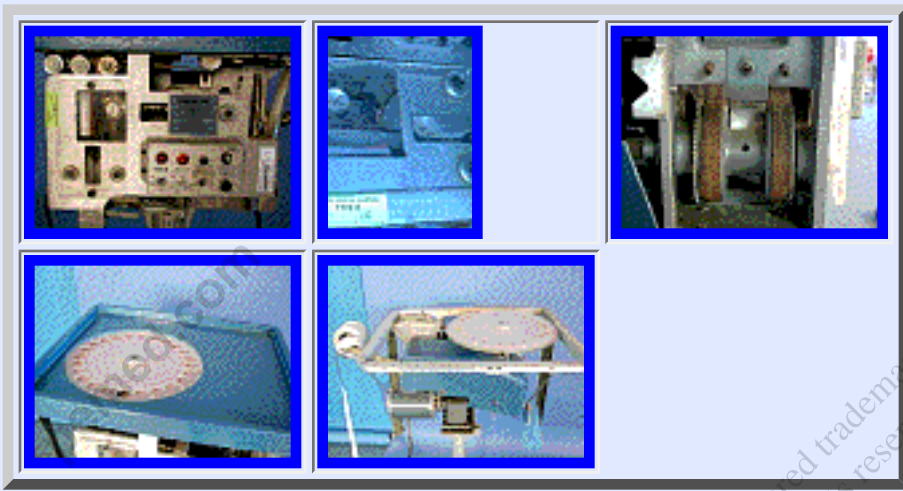

Shipping weight is about 200 lbs.


Bag / pouch sealer - \$675



List: \$3,200. MedCon: EX. Impulse bag sealer, made by Rennco, model "lift seal".
Seals bags up to 18" wide.

Unit will work with the following materials:

- Polyethylene bags
- Tyvek pouches
- paper pouches
- plastic pouches

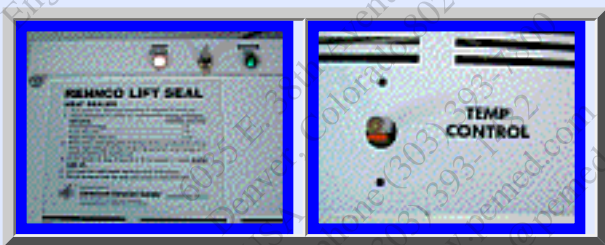

Exterior dimensions:

- length: 10"
- width: 19"
- height: 7"

Shipping weight is about 50 pounds each. 2 units in stock.

OmniCell Omnisupplier control module - \$1,500



List: \$45,600. MedCon: EX. Omnicell Omnisupplier control panel. This is a control module, the electronic guts of the system, for running a Omnicell optimized medication distribution system. May be used to dispense, control and automatically track supplies and medications in patient care areas. This unit is running V Jw4.1.1 software, 8/31/99, whatever that is. It is less than 1 year old, and appears new except has a crack in the plastic of the control panel.



For full information please see the manufacturers' web site, <http://www.omnicell.com/product/index2.asp>

Drug preparation cabinet - \$1,275



List: \$5,030. MedCon: G. Market Forge Drug preparation cabinet, model medi-prep. All stainless steel construction. Features:

- built in sink
- lap style locking drawer
- 3 shelf refrigerator

- under sinks service and storage area
 - 4 shelf, stepped vial display area
 - 2 adjustable height shelves
 - locking 3 shelf compartment
 - locking pull down, roll up overall cover.
 - dimensions: 7' high, 2' deep, 4' wide.
-

Impulse bag sealer - \$450



List: \$1,340. MedCon: G. Floor model bag sealer. Sure tight sealer, made by Electroven manufacturing. A floor model impulse type sealer, with foot control. Simple but reliable. Sealing bar is about 12" wide.

Shipping weight is about 100 pounds.

Liquid dispensers



- List: \$5,513. MedCon: G. Beckman Accuprep system 221 dispensing system, dispenses 2 vials at a time. - \$850
 - List: \$4,225. MedCon: F. Beckman Accuprep system 111, dispensing system. Dispenses liquid one vial at a time - \$450
-

Narcotics control cabinet / dispenser - \$1,800



List: \$9,000. MedCon: EX. Documed narcotics control cabinet. Locking cabinet for controlled substance distribution. Holds 14 rows of 20 columns, each of different drugs. Also has 7 internal storage drawers and cabinets for additional locked storage. Operates with a series of keys and levers so that exact pill counts are dispensed. Records on a paper tape what is dispensed

to whom.

© Copyright 1993 - 2011 by Production Engineering. PEMED® is a registered trademark of Production Engineering. MedCon is also a trademark of Production Engineering. All rights reserved. Contact PEMED for license requests. The usual copyright notices and disclaimers. Except for search engines or link listing, no part of this website, including derived information, may be republished without express written permission and payment of fees.

Production Engineering - Medical Equipment Division

6035 E. 38th Avenue, Denver, Colorado 80207 USA (303) 393-7800 FAX (303) 393-1482

6035 E. 38th Avenue
Denver, Colorado 80207
USA
telephone (303) 393-7800
FAX (303) 393-1482
<http://www.pemed.com>
email: markz@pemed.com

© Copyright 1993 - 2011 by Production Engineering. PEMED® is a registered trademark of Production Engineering. MedCon is also a trademark of Production Engineering. All rights reserved.

6035 E. 38th Avenue
Denver, Colorado 80207
USA
telephone (303) 393-7800
FAX (303) 393-1482
www.pemed.com
markz@pemed.com